


Terragene Disinfection Control - Thermal

Monitoring and checking the disinfection efficiency of moist heat/thermal disinfection in Washer Disinfectors

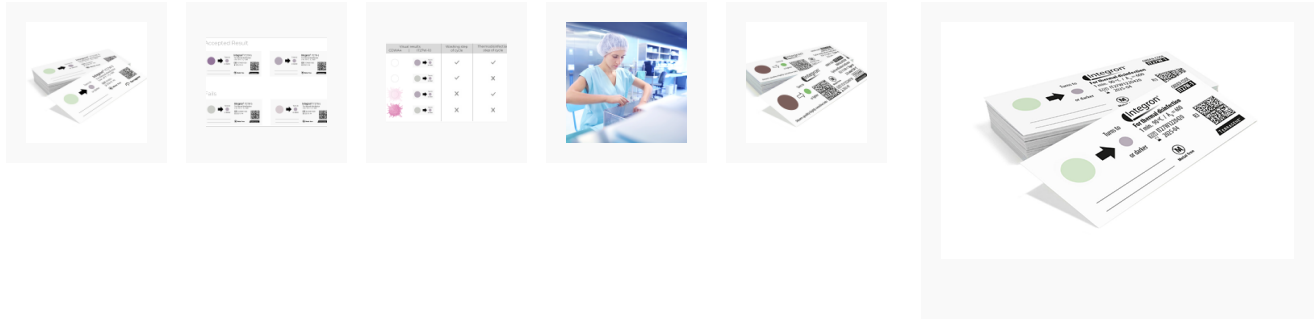

Features

- After cleaning, surgical instruments must be disinfected to ensure safe subsequent handling and processing.
- Thermal disinfection with moist heat is the most common method for disinfection of medical devices in the hospital setting.
- Terragene Integron Indicators have been designed to react in moist heat disinfection processes in washer disinfectors.
- Terragene Integron Indicators offer an accurate, convenient method of routine control, assuring disinfection efficiency.
- The green ink on the Terragene Integron Indicator has been developed to turn purple when the temperature and time conditions are met.
- Terragene Integrom IT27W-1 Indicator features a one side laminated strip, 70 x 30 mm, printed with indicator ink for routine control of thermal disinfection processes. Turns from Green to Purple.
- **1 minute, 90°C Program, A₀: 600** - Inactivates bacteria, mycobacteria, fungi and heat sensitive viruses (IT27W-1)
- **5 minutes, 90°C Program, A₀: 3000** - Efficacy against heat resistant viruses (Hepatitis B). Recommended by the Robert Koch Institute (IT27W-5)
- **10 minutes, 93°C Program, A₀ > 10.000** - Epidemic or “BGA/RKI” Program. High safety (IT27W-10)

Product Code	Product Description	Pack Size	Quote
IT27W-1	Terragene Integron Thermal Disinfection Indicator, 1 minute, 90°C, A ₀ : 600	200 Pk	
IT27W-5	Terragene Integron Thermal Disinfection Indicator, 5 minute, 90°C, A ₀ : 3000	200 Pk	
IT27W-10	Terragene Integron Thermal Disinfection Indicator, 10 minute, 93°C, A ₀ > 10.000	200 Pk	

Relevant Product

- [Smeg GW0260](#)
- [Smeg GW2145](#)
- [Smeg GW6010](#)
- [Smeg GW7010](#)

- [Smeg GW4060](#)
- [Smeg GW4190](#)
- [Smeg GW3060](#)
- [Smeg GW4090](#)
- [Smeg GW6290](#)
- [Smeg GW7015](#)
- [Belimed WD290](#)
- [Belimed PH810](#)
- [Belimed PH820 / PH840 / PH860](#)
- [Belimed PH880](#)

Applications

- Life Science
- Laboratory Infection Control
- CSD
- Healthcare Facilities
- Dental

Resources

- Video
- Video
- Brochure